

Micro-CT SSP

Measuring slice sensitivity profiles, slice thickness or axial modulation transfer function becomes easy with the Micro-CT-SSP-Phantom.

MicroCT-SSP contents a thin insert of material that highly attenuates X-ray (Au). The metal insert is only 10 microns thick and thus adapted to the nominal slice thickness of micro-CT scanners. The high-contrast platelet (diameter 1 mm) is embedded in the center of the cylindrical phantom made of resin.

Standard phantoms are scaled down from the diameter of an average human trunk (32 cm) to that of a fat mouse or lean rat (32 mm). Other dimensions as well as inserts with substances other than specified here can be produced upon request.

Specifications

Base material	resin
Diameter of gold insert	1 mm
Thickness of gold insert	10 μm

Dimensions for mouse-size phantom:

Phantom	diameter	32 mm
Phantom	length	. 40 mm
Phantom	weight	35 g


References:

[1] Kalender, W., Durkee, B., Langner, O., Stepina, E., Karolczak, M.: Comparative Evaluation: Acceptance Testing and Constancy Testing for Micro-CT Scanners.


Biomedizinische Technik 50 (2005), 1192-1193



Micro-CT SSP Phantom

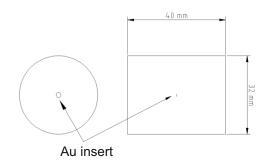

Schematic drawing of the phantom

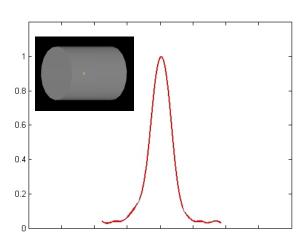

Example of scanner SSP

深圳市展业达鸿科技有限公司

吕先生:15920060912(微信同号) 0755-22934005(座机)

地址:深圳市福田区八卦二路八卦岭工业区615栋419 邮箱:hongqi@thingstet.com 网址:www.thingstest.com