


Fetal Ultrasound Biometrics Phantom

Model 068

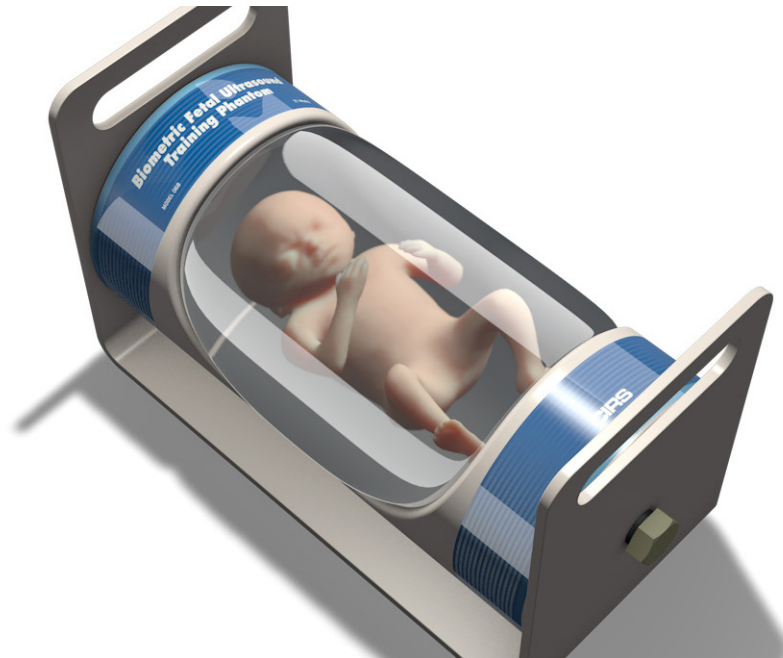

DEMONSTRATE • TEACH • PRACTICE

The mid-term ultrasound exam is used to estimate gestational age and to diagnose fetal malformation so students can gain competency with hands-on practice. Unfortunately, access to patients can be limited. The CIRS Model 068 Fetal Ultrasound Biometrics Phantom facilitates teaching and demonstration of fetal ultrasound techniques in a non-stressful situation.

A tissue equivalent full fetal model is suspended in a non-echoic, amniotic fluid-like environment. The model is housed in a rotatable cylinder with 2 fields of view. A variety of fetal/transducer orientations can be achieved.

The Model 068 contains an asymmetric head with upper portion of the skull, right and left brain lobes and lateral and third ventricles. These anatomical references are used to measure the biparietal diameter (BPD) and anterior/posterior diameter (APD). Right and left femoral shafts with distal epiphysis are provided for femur length (FL) measurements. An umbilical marker indicates proper position for taking abdominal circumference (AC). Crown-Rump length (CRL) can also be

taken. In addition, the model has full facial details making it suitable for demonstrating multi-dimensional systems and assessment of fetal anomalies.

All anatomies are based on published biometric data for normal fetal growth rates for a gestational age of 21 weeks. This enables assessment of composite measurement techniques and biometric analysis programs common to most ultrasound scanners.

Features

- Demonstrate system capabilities
- Teach imaging protocols to estimate gestational age
- Practice 2D, 3D & 4D fetal scan techniques
- Measure CRL, BPD, FL and AC
- Gain competency in performing prenatal ultrasound

深圳市展业达鸿科技有限公司

吕先生: 15920060912 (微信同号) 0755-22934005 (座机)

地址: 深圳市福田区八卦二路八卦岭工业区615栋419

邮箱: hongqi@thingstet.com

网址: www.thingstet.com



3D Reconstruction



Axial Brain & BPD Measurement

Sagittal Brain

Sagittal Abdomen with IVC Marker



Femoral Length for Measurements

Abdominal Circumference for Measurement

SPECIFICATIONS

DIMENSIONS:	35 cm x 17 cm x 17 cm (5.4" x 2.6" x 2.6")
TARGETED BIOMETRIC DIMENSIONS:	BPD: 4.2 cm (10th percentile) APD: 5.8 cm Head Circumference: 16 cm (16th percentile) Femoral Length: 3.8 cm Abdominal Circumference: 17.5 cm Crown-Rump Length: 22 cm
PHANTOM WEIGHT:	8.9 kg (24 lbs.)
MATERIALS:	Housing: PVC, Vinyl & ABS Phantom: Proprietary Gels

MODEL 068 INCLUDES

QTY	COMPONENT DESCRIPTION
1	Fetal Ultrasound Biometrics Phantom
1	Nylon Transducer Cover
1	User's Guide
-	12 Month Warranty*

*Please note: once any device has been inserted into the phantom (such as biopsy needles, localization wires, etc.) the 12 month warranty will not cover claims related to material desiccation or needle tracking.

OPTIONAL ACCESSORIES

MODEL	DESCRIPTION
9505	Case for Ultrasound Biometrics Phantom

©2013 Computerized Imaging Reference Systems Inc. All rights reserved.
Specifications subject to change without notice.
Publication: 068 DS 072120

Reference:

Hansmann, M. (1985) Ultrasound Diagnosis in Obstetrics and Gynecology, Springer, Berlin.

深圳市展业达鸿科技有限公司

吕先生: 15920060912 (微信同号) 0755-22934005 (座机)

Computerized Imaging Reference Systems, Inc. has been certified by UL DQS Inc. to ISO 13485:2016. Certificate Registration No. V50000005, MP2616

5栋419

邮箱: hongqi@thingstet.com

网址: www.thingstet.com